





Oral-B®

PROFESSIONAL
CARE

powered by **BRAUN**

500

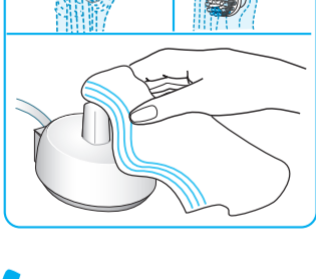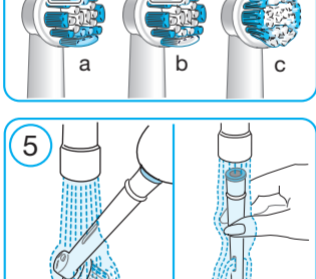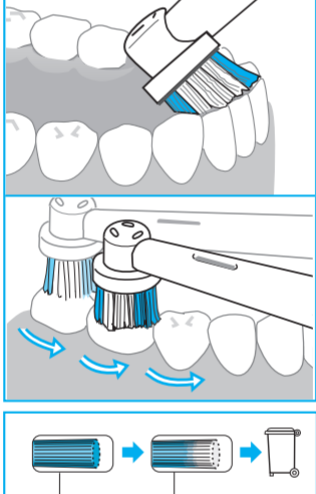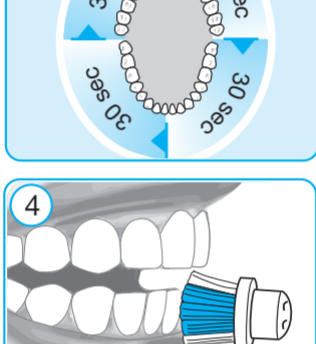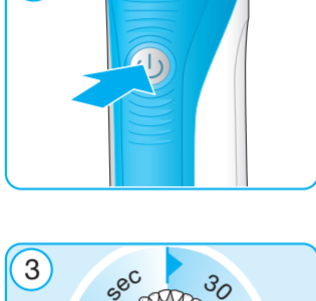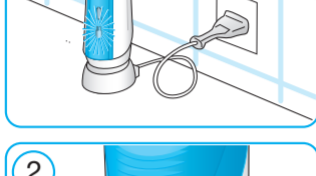- (GB)** 0 800 783 7010
- (IE)** 1 800 509 448
- (FR)** 0 800 944 803
(appel gratuit depuis un poste fixe)
- (BE)** 0 800 14 592
- (LU)** 800 21172
- (TR)** 0 800 261 63 65
- (RU)** 8-800-200-20-20
- (UA)** 0 800 505 000
- (ZA)** **Consumer Care Line**
0860 112 188
(Sharecall charged at local rates)

Internet:
www.oralb.com
www.braun.com
www.service.braun.com

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96199851/II-13
 GB/FR/TR/RU/UA/KZ/Arab

English

Thank you for choosing the Oral-B Professional Care 500. This toothbrush has been designed to meet the highest standards of quality, functionality and design.

Before operating this toothbrush, please read these instructions and save this manual for future reference.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for use by children under age three. Children from 3 to 14 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the toothbrush, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

- a Precision Clean brush head
- b On/off button
- c Handle
- d Charge indicator light
- e Low charge indicator light
- f Basic charger

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the basic charger (f) into an electrical outlet and place the toothbrush handle (c) on the charging unit.
- The green charge indicator light (d) flashes while the handle is being charged. Once it is fully charged, the green charge indicator light shines continuously for 5 seconds; then it turns off. A full charge enables up to seven days of regular brushing (twice a day, 2 minutes) (picture 1).
- When the rechargeable battery is low, the red low charge indicator light (e) starts flashing after turning off the toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Please note: If the charge indicator light is not flashing immediately, continue charging. It will appear in 10–15 minutes.

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Type: 3757

Noise level: ≤65 dB (A)

Using your toothbrush

Wet brush head (a) and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance (b). Slowly guide the brush head from tooth to tooth, spending a few seconds on each tooth surface. Brush both gums and teeth; first the outsides, then the insides and finally the chewing surfaces (picture 4). Brush at least 2 minutes to thoroughly remove plaque. Do not press too hard or scrub. Turn off the handle by pressing the on/off switch.

Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush all four quadrants of your mouth equally (picture 3). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorised when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets.


Pressure Sensor

To promote optimal brushing your toothbrush has a pressure control feature installed. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition you will also hear a different sound while brushing.

Individual brush heads and accessories

Oral-B Professional Care offers you a variety of brush heads and accessories which are available in stores that sell Oral-B products.

Oral-B Precision Clean brush head

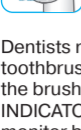cupping each tooth with a curved bristle trim and reaches deep between teeth at the same time.

Dentists recommend replacing any toothbrush every 3 months. Therefore the brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need (picture 4a). With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period, signaling it is time to change your brush head (picture 4b). If the bristles splay before the color recedes, you are exerting too much pressure on teeth and gums (picture 4c).

Cleaning recommendations

After brushing, rinse your brush head carefully under running water for some seconds with the handle switched on. Switch off the handle and remove the brush head. Clean both parts separately under running water; then wipe them dry before reassembling the toothbrush. Never place the charger in dishwasher or water; it should be cleaned with a damp cloth only (picture 5).

Subject to change without notice.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Service Centre.

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